



Accumax
Global Pty Ltd

Providing solutions to manage & reduce RISK!

Edition 3



MAXSafe

DUAL EYE

**PERSONAL
EYEWASH**

MAXBund

MAXSafe

MAXChock

DUAL EYE PERSONAL EYEWASH

WE COVER THE ENTIRE SPECTRUM

YOU ONLY HAVE TWO EYES – TAKE CARE OF THEM

Your eyesight is one of the most important senses. So take care of your eyes. Regardless how careful you are, accidents can happen, but everybody can make sure that they are prepared. With our eye wash solutions, you have done your best.

The green Saline Eye Wash deals with dirt. The blue pH Neutral is crucial in connection with accidents involving acids and alkali.

If an accident happens, seconds can decide whether your eyesight is damaged or not. Therefore, quick access to eye wash is crucial for all workplaces where there are risks of eye injuries.

MAXSafe 


A WIDE RANGE OF UNIQUE RINSING SOLUTIONS

Our wide product range ensures the best possible solution specifically for your needs. Choose small and large wall-mounted stations directly at the workplace. Or small, mobile solutions for the van, the tool box and the first aid box, and on the cleaning trolley in the special belt bag.



HOW TO PREVENT EYE INJURIES

Our concepts have been developed in close collaboration with skilled eye specialists. This means that the products have been thought out down to the smallest detail. They are easy to use with clear pictograms showing the correct use. The eye cup helps keep the eye open. The rinsing liquid hits the eye in the right place and with the right amount in a soft and even flow.

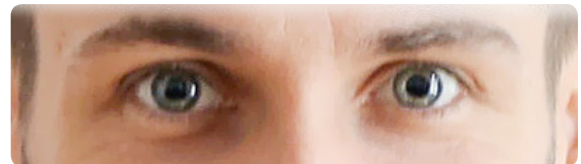


Saline Eye Wash removes foreign objects like dust and dirt from your eye, while pH Neutral will neutralize acids and alkali



Duo Eye Wash is easy to use if an accident happens



Turn the eye cup in the direction of the arrow until the seal breaks

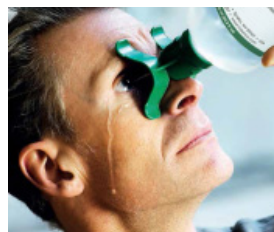

Lean back, squeeze the bottle and rinse



The eye can also be rinsed while you lean forward, by squeezing the bottle gently

1000ML SALINE EYE WASH DUO

Bottle with 0.9% sterile sodium chloride. Contains 1000 ml eye wash, which provides a longer rinsing effect. The ergonomic DUO eye cup makes it possible to rinse both eyes at the same time. The bottle can be used individually or as refill at DUO eye wash stations and boxes.

RINSING TIME: about 5 mins

PACKAGING: 1

PRODUCT CODE: AM-H7584



SALINE EYE WASH STATION DUO – WALL MOUNT

Wall station with 1000 ml Saline Eye Wash DUO. Particularly suitable for small and mobile workplaces where prolonged rinsing of both eyes may be required.

RINSING TIME:

about 5 mins

PACKAGING: 1

DIMENSIONS:

375mm H x 180mm W
x 95mm D

PRODUCT CODE: AM-H7543



500ML pH NEUTRAL EYE WASH DUO

A handy bottle with 4.9% sterile phosphate buffer that will quickly neutralize acids and alkali. The ergonomic DUO eye cup makes it possible to rinse both eyes at the same time. It is particularly suitable for workplaces where there is a risk that both eyes are affected. The bottle can be used individually or as refill at eye wash stations and boxes.

RINSING TIME: about 2 mins

PACKAGING: 1

PRODUCT CODE: AM-H7583



SALINE / pH COMBI-STATION DUO – WALL MOUNT

Wall station with 500 ml pH Neutral Eye Wash DUO and 1000 ml Saline Eye Wash DUO. Particularly suitable for places where there are risks for both eyes to be affected by foreign bodies, acids, and alkali.

RINSING TIME:

- pH Neutral Eye Wash Duo about 2 mins.
- Saline Eye Wash Duo about 5 mins

PACKAGING: 1

DIMENSIONS:

377mm H x 250mm W
x 95mm D

PRODUCT CODE:

AM-H7541



COMBI-BOX DUO MAXI – WALL MOUNT

Dust-proof wall station in impact-proof polystyrene with 500 ml pH Neutral Eye Wash DUO and 1000 ml Saline Eye Wash DUO. Particularly suitable for dusty production areas where there are risks for both eyes to be affected by accidents involving foreign bodies, acids and alkali.

RINSING TIME:

- pH Neutral Eye Wash Duo about 2 mins
- Saline Eye Wash Duo about 5 mins

PACKAGING: 1

DIMENSIONS: 357mm H x
252mm W x 116mm D

PRODUCT CODE: AM-H7554



ACCIDENTS INVOLVING ACIDS AND ALKALI WILL HAPPEN WHEN YOU LEAST EXPECT THEM



pH Neutral DUO saved the eyesight of an employee at a recycling plant. They handle plastic waste such as plastic bottles with cleaning detergents.

The liquid from a cleaning bottle went into both the employee's eyes.

Thanks to pH Neutral DUO, the harmful substance was neutralized, and the employee was back at work the next day.

EYE WASH IN CASE OF FOREIGN OBJECTS

Dual Eye Personal Eye Wash can remove discomfort caused by foreign bodies in your eyes and prevent subsequent complications.

Experience shows that about 50% of all eye accidents involving chemicals affect both eyes. With our DUO eye wash, you are prepared. You can rinse both eyes at the same time and will not have to choose which eye to save.

Time is a critical factor, especially with accidents involving acids/alkali. Of course the injury depends on the strength of the acid or alkali, but apart from that, the contact time is crucial. The optimal first aid must therefore neutralize the harmful substance as quickly as possible. With the unique DUO eye wash, you will treat both eyes at the same time.



About 50% of all accidents affect both eyes

We reserve the right to change the product specification without notice.

